



Europäisches Patentamt  
European Patent Office  
Office européen des brevets



(11) **EP 0 386 229 B2**

(12) **NEW EUROPEAN PATENT SPECIFICATION**

(45) Date of publication and mention  
of the opposition decision:  
**09.06.2004 Bulletin 2004/24**

(51) Int Cl.7: **C07H 21/00**  
**// C12Q1/68**

(45) Mention of the grant of the patent:  
**23.03.1994 Bulletin 1994/12**

(86) International application number:  
**PCT/GB1989/001114**

(21) Application number: **89910878.1**

(87) International publication number:  
**WO 1990/003382 (05.04.1990 Gazette 1990/08)**

(22) Date of filing: **21.09.1989**

(54) **SUPPORT-BOUND OLIGONUCLEOTIDES**  
SUBSTRAT-GEBUNDENE OLIGONUKLEOTIDE  
OLIGONUCLEOTIDES LIES A UN SUPPORT

(84) Designated Contracting States:  
**AT BE CH DE FR GB IT LI LU NL SE**

• **MASKOS, Uwe**  
**Oxford OX1 2JH (GB)**

(30) Priority: **21.09.1988 GB 8822228**

(74) Representative: **Hallybone, Huw George et al**  
**Carpmaels and Ransford,**  
**43 Bloomsbury Square**  
**London WC1A 2RA (GB)**

(43) Date of publication of application:  
**12.09.1990 Bulletin 1990/37**

(73) Proprietor: **OXFORD GENE TECHNOLOGY**  
**LIMITED**  
**Kidlington, Oxford OX5 2HB (GB)**

(56) References cited:  
**EP-A- 0 090 789**                    **EP-A- 0 241 363**  
**EP-A- 0 261 283**                    **EP-A- 0 305 929**  
**WO-A-85/01051**                    **DE-A- 3 446 635**

(72) Inventors:  
• **SOUTHERN, Edwin, M.**  
**Oxford OX1 3QU (GB)**

**EP 0 386 229 B2**

**Description**

**[0001]** There are several potential applications for oligonucleotides bound to solid supports. They could be used to test for the presence of mutations in complex DNAs - for example for disease loci in Humans. They could be used to select specific nucleic acids from the complex mixtures; for example specific mRNAs from a whole cell population. They will be useful in the invention described in International Application PCT/GB89/00460 filed 2 May 1989.

**[0002]** Several papers describe methods for attaching nucleic acids to solid matrices (1-6). These methods suffer from two problems: they require complex and often inefficient steps; the nucleotide is linked through the bases as well as the ends. Linkage through the bases interferes with subsequent use of the bound polynucleotide in hybridisation reactions.

**[0003]** Methods for synthesising oligonucleotides on solid supports are well established (7-14). The linkage between the oligonucleotide and the support is labile to the final reagent used to remove blocking groups in the bases, and so this step in the process also removes the oligonucleotide from the solid support. Oligonucleotides would remain tethered to the support if a stable link were used. Sproat and Brown (1985) have shown that a urethane link is more stable than the usual succinate link, however, it requires a complex synthesis, and is not completely stable to the final deprotection step.

**[0004]** Crea and Horn (1980) used a ribonucleotide, linked through the 5'-hydroxyl group to cellulose, to initiate oligonucleotide synthesis. The link between the first and second residues of the resulting chain is labile to the final deprotection step.

**[0005]** Arnold and Berg (1985) describe a polymeric support with a covalently bonded primer for oligonucleotide synthesis, wherein the primer is cleaved by selective oxidation without oxidizing other bonds of the oligonucleotide.

**[0006]** It is an object of this invention to provide a new link which is easy to synthesise and completely stable to standard deprotection steps.

**[0007]** In one aspect the invention provides a method of making a derivatised support suitable for oligonucleotide synthesis, which method comprises attaching a nucleoside 3'-phosphite reagent to a support carrying hydroxyl groups by a covalent phosphodiester link which is stable to conditions used for removing protective groups from oligonucleotide chains, characterised in that the hydroxyl groups are aliphatic hydroxyl groups in which the aliphatic moiety is hexaethoxy.

**[0008]** In another aspect the invention provides a method of preparing an oligonucleotide bound to a support by


- a) attaching a nucleoside 3'-phosphite reagent to a support,
- b) synthesising on the supported nucleoside an oligonucleotide chain including protecting groups, and
- c) removing the protecting groups from the oligonucleotide chain,

wherein the support carries hydroxyl groups, whereby in step a) the nucleoside becomes attached to the support by a covalent phosphodiester link which is stable to the conditions used in step c), characterised in that the hydroxyl groups are aliphatic hydroxyl groups in which the aliphatic moiety is hexaethoxy.

**[0009]** In a further aspect the invention provides a derivatised support suitable for oligonucleotide synthesis comprising a nucleoside linked through the 3'-position to a support by means of a covalent phosphodiester link of the structure -O-PY-O-, where Y is a protected or unprotected oxygen atom, characterised in that the link to the support is through an aliphatic moiety which is hexaethoxy.

**[0010]** The nature of the support is not critical to the invention. It may be massive or particulate and may for example be of derivatised silica gel or Kieselguhr-polydimethyl-acrylamide, or controlled-pore glass, or a plain glass surface. What is essential is that it carry aliphatic hydroxyl groups, and these in a form which are accessible for reaction with a nucleoside reagent. The hydroxyl groups may be part of a polymeric structure, which either constitutes the solid support or is derivatised onto a solid support.

**[0011]** The nature of the nucleoside reagent is not critical to the invention. Reagents commonly used in oligonucleotide synthesis may be used here. Preferably, the reagent is a phosphoramidite (7 and 8). The reagents and the product formed are indicated in the following reaction scheme.

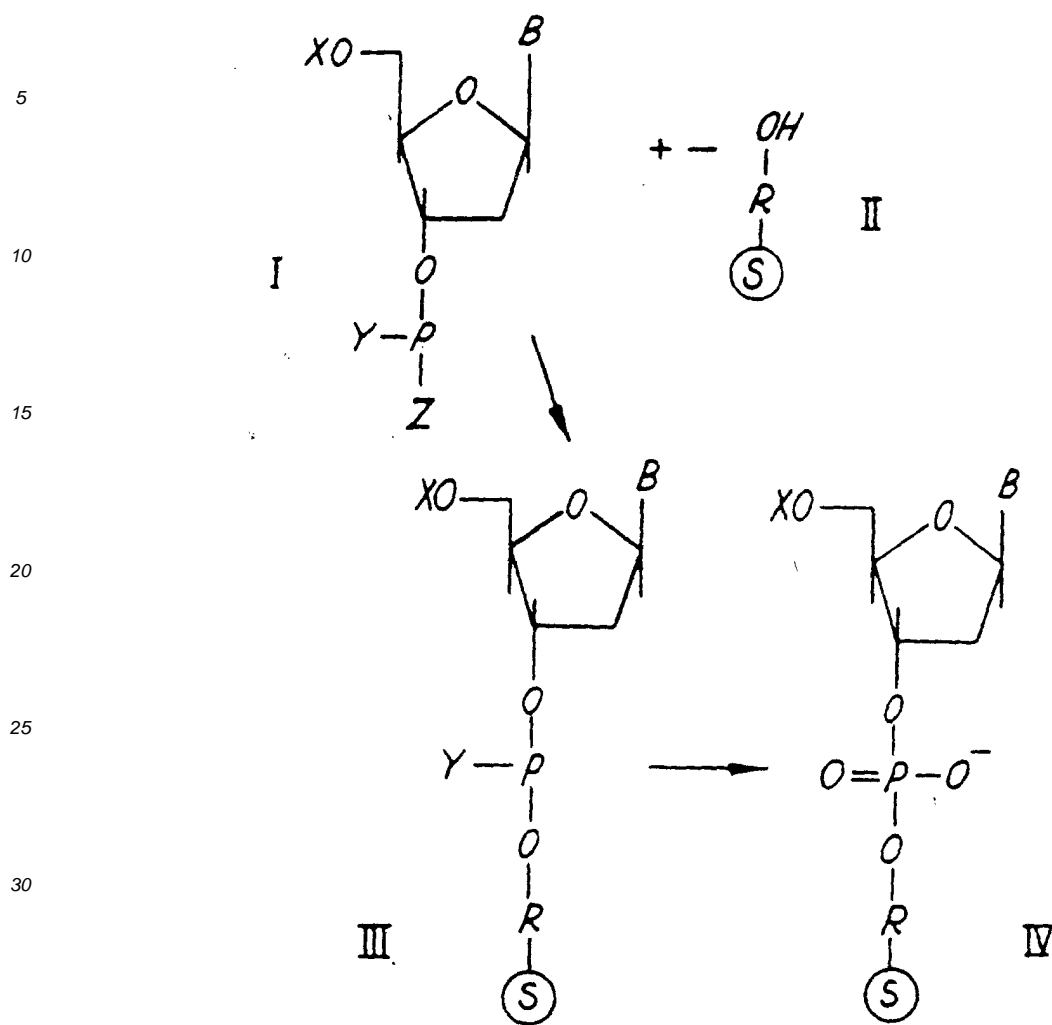

[0012] In this scheme, I represents the nucleoside 3'-phosphite reagent, II represents the derivatised support, III represents the covalent link initially formed by reaction between the two, and IV represents the final product.

[0013] B designates the base appropriate to the nucleoside concerned.

[0014] X may be a blocking group, such as a dimethoxytrityl group, whose nature is not critical to the invention; or may (particularly in IV) represent an oligonucleotide chain.

[0015] Y is a protected oxygen atom, generally an alkoxy group such as methoxy or beta-cyanoethoxy.

[0016] Z may be a di-(C1 to C4 alkyl)amine, or alternatively Cl or tetrazolyl.

[0017] R is aliphatic and may be an alkyl group.

[0018] S represents the solid support.

[0019] Automatic oligonucleotide synthesisers are commercially available. In prior techniques, it has been necessary to put into the synthesiser a solid support which has been pre-derivatised with the first nucleotide of the proposed oligonucleotide chain. This invention makes it possible to put into the synthesiser a solid support which carries aliphatic hydroxyl groups in which the aliphatic moiety is hexaethoxy. The nucleoside 3'- reagent is then attached to the support in the synthesiser as the first nucleotide of the proposed oligonucleotide chain.

[0020] In the above reaction schemes, conversion of III to IV involves an oxidation step. This may be effected using e.g. iodine or sulphur under standard conditions, either before or more usually after oligonucleotide synthesis.

[0021] The next step of the method involves synthesising on the supported nucleoside an oligonucleotide chain. Techniques for doing this are well known, and indeed automatic microprocessor-controlled machines are available to do the job. These techniques invariably involve the provision of protective groups to avoid unwanted side reactions, and a final step in the synthesis of any oligonucleotide involves removal of protective groups. It is this step that has previously resulted in solubilisation of the oligonucleotide. This step may typically involve removal of the methyl group from the phosphotriester groups e.g. by using thiophenoxide; removal of protecting groups from N atoms of the nucle-

otide bases, e.g. by means of ammonia at elevated temperature; and removal of a protective group from the 5'-position of the last nucleotide to have been added to the oligonucleotide chain, e.g. by means of dichloro acetic acid. The described covalent link between the initial nucleoside and the support is stable to these and other conventional deprotection steps.

5 **[0022]** The following examples illustrate the invention.

#### Example 1

10 **[0023]** Ballotini glass beads (20g, 90-130  $\mu\text{m}$  diam., Jencons), were suspended in a mixture of xylene (40ml), glycidoxypropyltrimethoxysilane, and a trace of diisopropylethylamine at 90°C overnight with stirring, then washed thoroughly with methanol, ether and air-dried. These derivatised beads (6g) were heated with stirring in hexaethyleneglycol containing a catalytic amount of concentrated sulphuric acid, overnight in an atmosphere of argon, at 80°C, to yield alkyl hydroxyl derivatised beads. After washing with methanol and ether, the beads were dried under vacuum and stored under argon at -20°C.

15 **[0024]** A small amount of the hydroxyalkyl derivatised beads was put into the reaction vessel of an automatic oligonucleotide synthesiser (Applied Biosystems), programmed to synthesise the sequence of the left cohesive end of bacteriophage lambda. The first nucleotide was a 3'-phosphoramidite (I) in which Y was beta-cyanoethoxy and Z was di-isopropylamine. This became covalently attached to the derivatised beads.

20 **[0025]** Each step in the synthesis can be monitored by measuring the amounts of trityl group removed, in the spectrophotometer. By this test, the stepwise yield was 96-99%. Thus both yield and purity were high, and we calculate that the 140mg of beads holds more than 4 nmol of the oligonucleotide.

**[0026]** The product was deprotected by the standard treatment with hot ammonia and washed thoroughly with distilled water. It was then used in a hybridisation with the complementary oligonucleotide cosL which had been labelled at the 5' end with  $^{32}\text{P}$ .

25 **[0027]** To 1mg of beads derivatised with the left end of phage lambda (cosL) was added the complementary oligonucleotide, radioactively labelled (13,000 cpm  $^{32}\text{P}$  in 20  $\mu\text{l}$ , 100 mM NaCl, 1mM EDTA). The mixture was incubated for 2 hours at 30 °C, the radioactive solution removed and the beads washed thoroughly with ice cold buffer. The oligonucleotide which remained attached to the beads (ca. 3,000 cpm, 23%) could be quantitatively removed by elution with distilled water.

30 **[0028]** A control "hybridisation" under identical conditions, with a non-complementary oligonucleotide showed 0.02% of non-specific binding to the beads.

**[0029]** These experiments show that the synthesis is straightforward and that the beads can be used successfully in hybridisation tests.

**[0030]** The glass beads are ideal for packing columns to provide an affinity matrix with many desirable properties.

35 **[0031]** The following Table 1 summarises physical parameters of beads derivatised by the general technique described in Example 1.

TABLE 1

Properties of derivatised ballotini beads		
Bead size (microns)	90 - 130	3 - 6
surface area per bead ( $\text{mm}^2$ )	0.045	0.00011
volume per bead ( $\text{mm}^3$ )	$9 \times 10^{-4}$	$1.1 \times 10^{-7}$
density ( $\text{g}/\text{cm}^3$ )	2.5	2.5
45 mass per bead (mg)	0.0023	$2.8 \times 10^{-7}$
number of beads per (mg)	435	$3.5 \times 10^6$
surface area per mg ( $\text{mm}^2$ )	19.5	385
oligonucleotide loading per mg (pmol)	~80	~80
50 oligonucleotide loading per bead (fmol)	184	0.02
values calculated for radius (microns)	60	3

#### Example 2

55 **[0032]** Two glass plates were clamped together with a narrow gap between them. Derivatisation was effected by means of the same reagents, which were injected into the narrow gap, and under the same conditions as Example 1.

**[0033]** Oligonucleotide synthesis was performed by hand under standard conditions using the derivatised glass plate

as a solid support. The first nucleotide was a 3'- hydrogen phosphate, used in the form of the triethylammonium salt. The yield and purity of both the first and subsequent steps of the oligonucleotide synthesis were high.

**[0034]** The following Examples demonstrate several ways of hybridizing labelled DNA fragments and oligonucleotides to the derivatised beads.

5

### Example 3

**[0035]** Hybridisation to glass beads attached to sticks.

10

**[0036]** A plastic stick 2cm long was dipped into molten polypropylene and then brought into contact with a pile of derivatised glass beads and allowed to cool. Approximately 100 - 200 beads adhered to the stick. This method of holding the beads greatly facilitates hybridisation as will be shown in a number of typical experiments:

**[0037]** A stick with approximately 100 glass beads derivatised with the sequence 3' AGG TCG CCG CCC 5' was dipped into 30  $\mu$ l of a solution containing 0.1 M NaCl and 80 fmol of the complementary oligonucleotide, labelled at the 5' end to an activity of 30,000 cpm.

15

**[0038]** After 30 min at 30 °C the stick was removed from the tube, rinsed and the bound material eluted by dipping the stick into 0.1 M NaCl at 50 °C. The amount of oligonucleotide hybridised was then determined by scintillation counting.

**[0039]** Typically 4% of the input oligonucleotide could be picked up this way; 0.1% were bound nonspecifically and could not be removed. Thus the binding capacity of a single bead is approximately 0.03 fmol oligonucleotide.

20

**[0040]** Larger proportions of the initial oligonucleotide could be picked up by decreasing the temperature and increasing the length of the hybridisation. Thus in a similar experiment 5.3% was hybridised at 30 °C after 55 minutes with 0.05% bound nonspecifically, and 13% after 16 hours at 30°C with 0.2% binding nonspecifically.

**[0041]** As a control, noncomplementary oligonucleotide 5' GGG CGG CGA CCT 3' showed only 0.2% binding after 14 hours.

25

**[0042]** In summary, experiments with derivatised glass beads attached to plastic sticks have proved to be very easy and shown the high specificity of hybridisation to the beads.

### Example 4

30

Batch hybridisation.

**[0043]** The amount of radioactive material hybridising to the beads could be increased further, and the non-specific binding decreased by carrying out hybridisation with beads typically 1 mg in 0.5 ml centrifuge tubes. After hybridisation the tubes were spun, the supernatant removed and the beads washed. Washing at a temperature higher than  $T_m$  resulted in complete melting of the hybrids so that the bound material could be measured by Cerenkov counting. In this way we determined the dependence of rate of hybridisation and elution on salt concentration and temperature as follows:

35

**[0044]** To each of 5 tubes is added an approximately equal number of beads and complementary oligonucleotide (30,000 cpm) in 50  $\mu$ l of 0.1M NaCl. Hybridisation was carried out at 30°C overnight to maximise the amount of hybrid. The solution was removed and the beads washed twice with ice-cold 0.1 M NaCl. Elution was for increasingly longer times at a different temperature for each tube. After each interval the supernatant was removed, the beads were washed twice with ice-cold NaCl solution (100  $\mu$ l 0.1 M), and eluted with prewarmed NaCl solution (100  $\mu$ l 0.1M).

40

**[0045]** Table 2 details the percentage of bound oligonucleotide eluted in the course of time. There is a clear dependence of elution rate on temperature. For example, three times as much material eluted at 65 °C than at 30 °C, within 5 minutes. Not surprisingly, even at 30 °C which is well below  $T_m$ , there is a non-negligible rate.

45

**[0046]** In another experiment the concentration of input oligonucleotide was varied 50 fold. Hybridisation for 2 hours at 30°C in 0.1M NaCl, 1mM EDTA in 1.5ml centrifuge tubes was followed by removal of the supernatant, three washes with 0.1M NaCl at 0°C and the amount of radioactivity associated with the beads determined by Cerenkov counting. There is an almost linear relationship between concentration and amount of hybrid (i.e. rate of hybridisation, Table 3), which suggests that the hybridisation is a pseudo first order reaction.

50

**[0047]** The highest concentration of oligonucleotide was 160 fmol (corresponding to 13,000 cpm) in 20  $\mu$ l. Only 0.03% bound nonspecifically and could not be eluted.

**[0048]** Furthermore, only 0.02% non-complementary oligonucleotide bound to the beads in a similar experiment (12 out of 65,000 cpm), a further indication that this method of isolating DNA fragments is very specific and clean.

55

TABLE 2

Effect of temperature on elution rate						
Elution time temperature	5'	20'	50'	110'	330'	% remaining
	% eluted					
30°C	23	38	54	71	85	0
36°C	34	61	78	90	96	3
44°C	42	54	67	87	97	2
65°C	74	89	96	98	98	1.4

TABLE 3

Effect of oligonucleotide concentration on hybridisation rate					
Relative concentration	50	20	10	2	1
% hybridised	23	14	5.5	0.8	0.6
relative rate	38	23	9	1.5	1

Example 5

Isolation of longer DNA fragments:

**[0049]** The isolation of longer DNA fragments was demonstrated in the following experiments: Total  $\lambda$  DNA (5  $\mu$ g in 35  $\mu$ l restriction enzyme buffer) was digested with 15 units of the restriction endonuclease Hinf 1 and the resulting 143 restriction fragments dephosphorylated with calf intestinal phosphatase (1 unit). After one hour at 37°C 4  $\mu$ l of EGTA was added, the mixture incubated for another 45 minutes at 65 °C, phenol extracted and ethanol precipitated.

**[0050]** The mixture was dissolved in 50  $\mu$ l kinase-labelling buffer (8  $\mu$ l 10 x PNK buffer, 1  $\mu$ l 0.1 M DTT, 4  $\mu$ l PNK enzyme, 30  $\mu$ l distilled water, 15  $\mu$ Ci  $^{-32}$ P.ATP), incubated for one hour at 37°C, made up to 100  $\mu$ l and spun down a Sephadex G25 column to remove non-incorporated nucleotides.

**[0051]** An aliquot of the ballotini glass beads (100  $\mu$ m, Jencons) was derivatised with the sequence of the right cohesive end of bacteriophage lambda, viz. 3' - CCC GCC GCT GGA 5', deprotected in ammonia, washed, and a small amount put at the bottom of a U-shaped capillary. This bottom part was kept at 40°C in a controllable temperature block, the upper left and right arms of the capillary were kept at 67°C by jacketing them with a plastic syringe and pumping hot water through it. The hybridisation solution was added (75  $\mu$ l 0.1 M NaCl containing 10 pmol 5' ends, 33 fmol left  $\lambda$  cohesive ends complementary to the oligonucleotide on the beads, total radioactivity ca. 700,000 cpm).

**[0052]** A pump was attached to one arm of the capillary and the hybridising solution cycled back and forth between the parts of the capillary that were kept at 40°C and 67°C respectively. Hybridisation of the left end to the beads would occur at 40°C and at 67°C the two sticky lambda ends that reannealed in solution would be denatured.

**[0053]** After 4 hours the hybridisation solution was removed, the beads washed extensively, then eluted in hot TE. An aliquot of the solution was loaded onto a 5% polyacrylamide gel. Autoradiography revealed only one band of the correct size in the washes and elution lanes. Altogether 1000 cpm (ca. 1.5 fmol) of the left end were bound by the beads which corresponded to 5% of the theoretical amount.

**[0054]** In summary, this experiment demonstrates the highly specific isolation of a long DNA fragment from a complex mixture.

Example 6

Column Chromatography.

**[0055]** Another easy and convenient way to isolate oligonucleotide and to test the hybridisation behaviour of the novel support is by column chromatography.

**[0056]** A glass capillary (diameter 1.0 mm) was drawn out at one end so as to yield a very narrow pointed opening. This was then plugged by filling in crushed glass particles from the other end and sintering them in the flame of a Bunsen burner so that they adhered to the glass. The inside of the capillary was silanised by passing through a solution of dichlorodimethyl-silane in trichloroethane and washing with ethanol. Approximately 40 mg of the glass beads were

layered on the glass frit and the top of the column connected to a syringe that could be driven by an infusion pump. In this way radioactive hybridisation solution and washing solutions could be applied to the column at different rates, typically in the range of 3 - 10  $\mu\text{l}/\text{min}$ .

**[0057]** In a typical experiment, 0.2 pmol of labelled oligonucleotide (34,000 cpm) in 1 ml of a solution of 0.1M NaCl, 0.1% SDS in TE pH 7.5 was applied to the column at a rate of 3  $\mu\text{l}/\text{min}$ . The jacketed column was kept at 35°C. 90  $\mu\text{l}$  fractions were collected in micro-centrifuge tubes and the amount of radioactivity determined by Cerenkov counting. A 0.1M NaCl washing solution was applied in the same way and collected. Raising the temperature in the jacket allowed us to recover the oligonucleotide.

**[0058]** Thus it was determined that 70% of the oligonucleotide bound to the glass beads and could be eluted at a higher temperature with only 0.1% of the material remaining on the support.

**[0059]** A control experiment with non-complementary oligonucleotide (mismatch at position 7) showed a remaining 400 cpm out of 80,000 applied after washing at 35°C. At 40°C only 140 cpm = 0.2% remained.

**[0060]** The percentage of accessible oligonucleotide on the support was determined. 0.13 pmol kinase-labelled and 5 pmol unlabelled oligonucleotide were applied to 40 mg beads. 10% of the material (ca. 500 fmol) hybridised to the support that contained a total of ca.3 nmol oligonucleotide, measured from the de-trityylation during synthesis.

**[0061]** From this result we calculate that one in 6000 oligonucleotides on the support hybridised under conditions used (35°C, low salt) where negligible binding of mismatched oligonucleotides occurs. The melting was very sharp again, with most of the oligonucleotides eluted in two 90  $\mu\text{l}$  fractions at 48 °C.

**[0062]** These experiments suggest that the derivatised beads will be useful in the chromatographic separation of nucleic acids.

## REFERENCES

### [0063]

1. Langdale J A and Malcolm A D (1985)  
A rapid method of gene detection using DNA bound to Sephacryl. *Gene* 1985, 36(3), 201-210.
2. Seed, B (1982)  
Diazotizable aniline cellulose paper for the coupling and hybridization of nucleic acids. *Nucl. Acids Res.* 10, 1799 - 1810.
3. Allfrey and Inoue (1978)  
Affinity chromatography of DNA-binding proteins on DNA covalently attached to solid supports. *Methods Cell Biol.* 17, 253 - 270.
4. Banemann, H; Westhoff, P and Herrmann G (1982)  
Immobilisation of Denatured DNA to Macroporous Supports: 1 Efficiency of different coupling procedures. *Nucl. Acids Res.* 10, 7163-7180.
5. Astell C R and Smith M (1972)  
Synthesis and Properties of Oligonucleotide-Cellulose Columns. *Biochemistry* 11, 4114-4120.
6. Gilham P.T., *Biochemistry*, 7, No 8, 1968, 2809-13.
7. Matteucci M D and Caruthers M H  
*J.Am. Chem.Soc.* 1981, 103, 3185-3191.
8. Beaucage S L and Caruthers M H  
*Tetrahedron Letters*, Vol 22, No.20, pp 1859-1862, 1981.
9. Adams S P et al, *J.Am. Chem.Soc.* 1983, 105, 661-663
10. Sproat D S and Brown D M  
*Nucleic Acids Research*, Vol 13, No.8, 1985, 2979-2987.
11. Crea R. and Horn T, *Nucleic Acids Research*, 8, No 10, 1980, 2331-48.
12. Andrus A. et al., *Tetrahedron letters*, Vol 29, No. 8, pp 861-4, 1988.
13. Applied Biosystems User Bulletin, Issue No. 43, Oct 1 1987, "Methyl phosphonamidite reagents and the synthesis and purification of methyl phosphonate analogs of DNA".
14. Miller P.S. et al., *Nucleic Acids Research*, 11, pages 6225-6242, 1983.
15. Arnold L.J. and Berg R.P., *Molecular Biosystems Inc PCT Application WO 85/01051*.

## Claims

1. A method of making a derivatised support suitable for oligonucleotide synthesis, which method comprises attaching a nucleoside 3'-phosphite reagent to a support carrying hydroxyl groups by a covalent phosphodiester link which is stable to conditions used for removing protective groups from oligonucleotide chains, **characterised in that** the hydroxyl groups are aliphatic hydroxyl groups in which the aliphatic moiety is hexaethoxy.

2. A method of preparing an oligonucleotide bound to a support by

- a) attaching a nucleoside 3'-phosphite reagent to a support,
- b) synthesising on the supported nucleoside an oligonucleotide chain including protecting groups, and
- c) removing the protecting groups from the oligonucleotide chain,

wherein the support carries hydroxyl groups, whereby in step a) the nucleoside becomes attached to the support by a covalent phosphodiester link which is stable to the conditions used in step c), **characterised in that** the hydroxyl groups are aliphatic hydroxyl groups in which the aliphatic moiety is hexaethoxy.

3. A method as claimed in claim 1 or claim 2, wherein the support is of porous glass or glass beads.

4. A method as claimed in any one of claims 1 to 3, wherein the nucleoside reagent is a nucleoside 3'-phosphoramidite.

5. A method as claimed in any one of claims 2 to 4, wherein step c) involves removing protecting groups from phosphotriester groups and N atoms of the nucleotide bases and from the 5'-position on the last nucleotide of the chain.

6. A derivatised support suitable for oligonucleotide synthesis comprising a nucleoside linked through the 3'-position to a support by means of a covalent phosphodiester link of the structure -O-PY-O-, where Y is a protected or unprotected oxygen atom, **characterised in that** the link to the support is through an aliphatic moiety which is hexaethoxy.

7. A derivatised support as claimed in claim 6, wherein the support is of porous glass or glass beads.

## Patentansprüche

1. Verfahren zur Herstellung eines derivatisierten Trägers, geeignet zur Oligonukleotidsynthese, wobei das Verfahren umfaßt, Binden eines Nucleosid-3'-phosphitreagenzes an einen Hydroxylgruppen-tragenden Träger durch eine kovalente Phosphodiesterbindung, die unter den zur Entfernung der Schutzgruppen für Oligonukleotidketten verwendeten Bedingungen stabil ist, **dadurch gekennzeichnet, daß** die Hydroxylgruppen aliphatische Hydroxylgruppen sind, in welchen der aliphatische Teil Hexaethoxy ist.

2. Verfahren zur Herstellung eines an einen Träger gebundenen Oligonukleotids durch

- a) Binden eines Nucleosid-3'-phosphitreagenzes an einen Träger,
- b) Synthetisieren einer Oligonukleotidkette, die Schutzgruppen einschließt, an dem befestigten Nucleosid und
- c) Entfernen der Schutzgruppen aus der Oligonukleotidkette,

wobei der Träger Hydroxylgruppen trägt, wodurch in Schritt a) das Nucleosid an den Träger über eine kovalente Phosphodiesterbindung gebunden wird, die unter den in Schritt c) angewendeten Bedingungen stabil ist, **dadurch gekennzeichnet, daß** die Hydroxylgruppen aliphatische Hydroxylgruppen sind, in welchen der aliphatische Teil Hexaethoxy ist.

3. Verfahren nach Anspruch 1 oder Anspruch 2, wobei der Träger poröses Glas oder poröse Glaskugeln darstellt.

4. Verfahren nach einem der Ansprüche 1 bis 3, wobei das Nucleosidreagenz ein Nucleosid-3'-phosphoramidit ist.

5. Verfahren nach einem der Ansprüche 2 bis 4, wobei Schritt c) die Entfernung von Schutzgruppen aus Phosphotriestergruppen und N-Atomen der Nucleotidbasen und aus der 5'-Stellung am letzten Nucleotid der Kette ein-

schließt.

- 5 6. Derivatisierter Träger, geeignet zur Oligonukleotidsynthese, umfassend ein Nucleosid, gebunden über die 3'-Stellung an einen Träger mit Hilfe einer kovalenten Phosphodiesterbindung der Struktur -O-PY-O-, wobei Y ein geschütztes oder ungeschütztes Sauerstoffatom darstellt, **dadurch gekennzeichnet, daß** die Bindung an den Träger über einen aliphatischen Teil, welcher Hexaethoxy ist, erfolgt.
7. Derivatisierter Träger nach Anspruch 6, wobei der Träger poröses Glas oder Glaskugeln darstellt.

10 **Revendications**

- 15 1. Procédé de fabrication d'un support fonctionnalisé convenant pour la synthèse d'oligonucléotides, lequel procédé comprend la fixation d'un réactif de type nucléoside 3'-phosphite sur un support portant des groupes hydroxy par une liaison phosphodiester covalente qui est stable dans les conditions utilisées pour l'élimination de groupes protecteurs de chaînes oligonucléotidiques, **caractérisé en ce que** les groupes hydroxy sont des groupes hydroxy aliphatiques dans lesquels le fragment aliphatique est un groupe hexaéthoxy.
- 20 2. Procédé de préparation d'un oligonucléotide fixé à un support, par
- a) fixation d'un réactif de type nucléoside 3'-phosphite sur un support,
- b) synthèse sur le nucléoside fixé sur le support, d'une chaîne oligonucléotidique comportant des groupes protecteurs, et
- 25 c) élimination des groupes protecteurs de la chaîne oligonucléotidique, le support portant des groupes hydroxy, de sorte que dans l'étape a) le nucléoside se fixe sur le support par une liaison phosphodiester covalente qui est stable dans les conditions utilisées dans l'étape c), **caractérisé en ce que** les groupes hydroxy sont des groupes hydroxy aliphatiques dans lesquels le fragment aliphatique est un groupe hexaéthoxy.
- 30 3. Procédé selon la revendication 1 ou 2, dans lequel le support consiste en verre poreux ou en perles de verre.
4. Procédé selon l'une quelconque des revendications 1 à 3, dans lequel le réactif nucléosidique est un nucléoside 3'-phosphoramidite.
- 35 5. Procédé selon l'une quelconque des revendications 2 à 4, dans lequel l'étape c) comporte une élimination de groupes protecteurs de groupes phosphotriester et d'atomes d'azote des bases nucléotidiques et de la position 5' sur le dernier nucléotide de la chaîne.
- 40 6. Support fonctionnalisé convenant pour la synthèse oligonucléotides, comprenant un nucléoside lié par la position 3' à un support au moyen d'une liaison phosphodiester covalente de formule -O-PY-O-, dans laquelle Y est un atome d'oxygène protégé ou non protégé, **caractérisé en ce que** la liaison au support est par un fragment aliphatique qui est un groupe hexaéthoxy.
- 45 7. Support fonctionnalisé selon la revendication 6, dans lequel le support consiste en verre poreux ou en perles de verre.